

1000X 1000 1000 1000 LCD 1000 1000 4.3 1000 1000 LCD 1000 1000

# Digital Electron Microscope

2 million pixel lens

50-1000X magnification

Connect to computer



# Digital Electron Microscope

2 MP | 50-1000X magnification | Support PC



4.3 inch  
screen



2 Million  
Pixels



50X-1000X  
magnification



Connect to PC

# 4.3-inch HD IPS screen

HD screen, exquisite picture quality



4.3 inch  
screen



178° viewing  
angle



1080P  
resolution



# Ultra-precision observation/imaging

2 million HD lens, 50-1000X magnification





50X



500X



1000X






# Stable storage and fast performance


The device has its own storage function, so you can take photos and record them without a memory card.



Photo

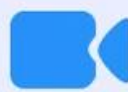

Video

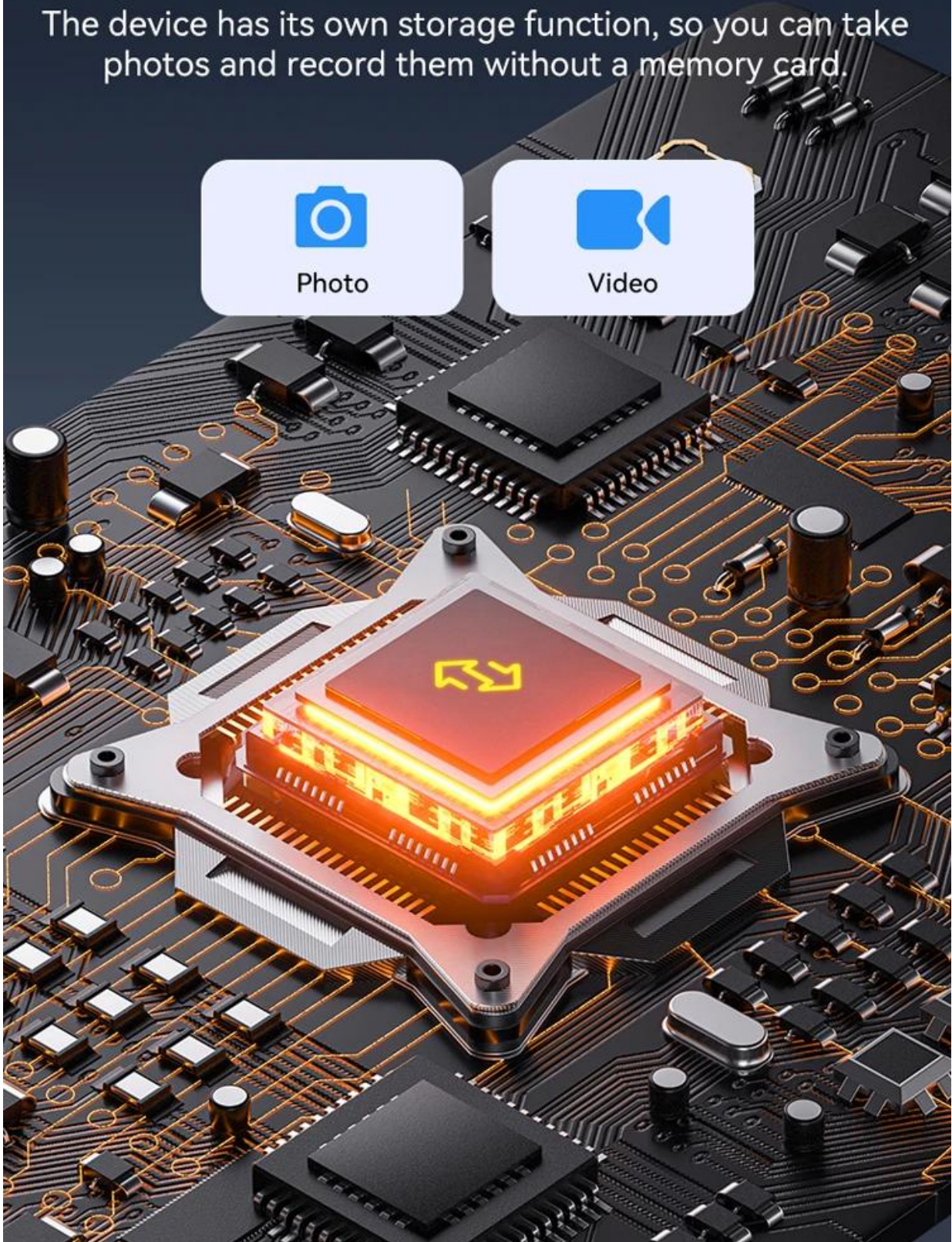



# Infinitely adjustable ring light

8 high-brightness LED ring lights illuminate  
dark environments



BRIGHTNESS ADJUSTMENT 0-100%:



# Can be connected to computer/laptop

Real-time transmission of HD images without delay



Mac



Windows





# Multi-scenario application

Microscope with screen, not limited to personal observation



Experimental observations



Classroom observation



Share screen




Observe flora and fauna

# Support taking photos and videos

Supports taking photos and videos, and can read TF card at any time for photo or video sharing



Support camera function

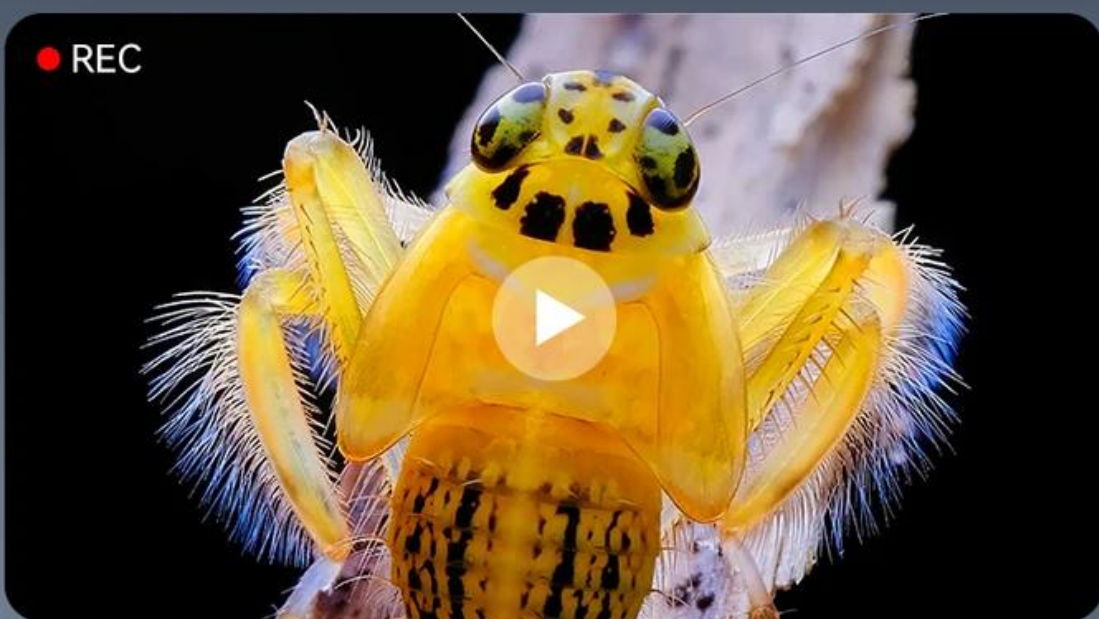

Support video recording function



# Product details

## TYPE-C transmission

Can be used to connect to a computer to transfer images and charge quickly



## Liftable storage table

Adjustment knob to adjust the vertical position of the shelf





## Silicone non-slip base


Stands stably on the desktop and is not easy to shake





# Operation guide


- 

1. Insert the MicroTF card (not equipped with this unit) into the MicroTF card slot; press gently until the card is locked ( Note: MicroTF cards must be inserted to record and store pictures and videos ).
- 

2. Long press on/off button, turn on the device.
- 

3. Put object to be observed onto the object stage, Rotate Wheel A to bring the object near to the microscope tube.
- 

4. Rotate the wheel B left and right to obtain fine focus.
- 

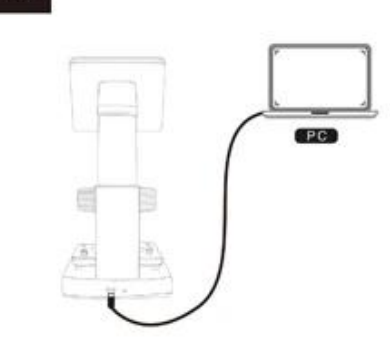5. When a clear focus got, rotate Wheel C clockwise to lock the object stage.
- 

6. Connect the microscope to your computer via a USB cable, turn on the device, select the "mass storage" mode, and view the stored photos and videos on your computer.





- |                                     |                          |                              |
|-------------------------------------|--------------------------|------------------------------|
| <b>1</b> 4.3*LCD display            | <b>6</b> Microscope tube | <b>11</b> Reset hole         |
| <b>2</b> MicroTF card slot          | <b>7</b> Wheel C         | <b>12</b> DC 5V in           |
| <b>3</b> USB power output port      | <b>8</b> Wheel A         | <b>13</b> Charging indicator |
| <b>4</b> LED brightness dimmer      | <b>9</b> Object stage    |                              |
| <b>5</b> Focusing wheel ( Wheel B ) | <b>10</b> Capture button |                              |

## Product size



□□

□□

□□ □□

□□

□□□□

□□

□□

□□

□□

□□

□□ □□

□□ □□

□□□ □□

□□□ □□

□

SE-I306

□□□ □□□

□□

□□□□

30X19.5X13.8cm

2MP

1920\*1080

1000X

LED 8□

0-40mm

JPEG

MP4/AVI



□□□ □□

□□ □□□

□□

□□ □ □□

1200mAh/1.5□□

WindowsXP, □□□, Win7/8/10, Mac10.5

CE/RoHS/ISO9001

# Product list/Parameters

Product packaging\*1, Device\*1, Instruction manual\*1,  
Cleaning cloth\*1, TYPE-C data cable\*1



Product name:  
**Digital microscope**

Product number:  
**306-B**

Screen:  
**4.3-inch IPS**

Pixels:  
**2 million**

Resolution:  
**1920\*1080P**

Gain:  
**50-1000X**

Best focal length:  
**0-40mm**

Light:  
**8 LEDs**

Port:  
**TYPE-C**

Battery capacity:  
**1200mAh**

Battery life:  
**1.5 hours**

Charging time:  
**2 hours**

## COMPANY INFORMATION

• SELECTECH •

About us



✓ Quality



✓ Experience



✓ Service



✓ Produce



### Selectech Electronics Co.,Ltd

- Welcome to ShenZhen Selectech Electronics Co., Limited. -- Your power sourcing for smart home and home health&beauty care products, consumer electronics, plastics products! Selectech is 11 years gold supplier in alibaba with trade assurance.



Office room



Team picture (1)



Team picture (2)

- Selectech was founded in Shenzhen in year 2005. Selectech invested factories focuses on production and R&D for smart home, home health care Otoscope&Camera,home consumer electronics. like Alexa Echo Google home control WiFi smart socket plug, USB charger wall socket receptable,smart fast wall charger, universal travel adapter,ENT otoscope,oral camera, USB&Wireless microscope and endoscope,DLP projector,baby monitor ect.



## Our factory



- Welcome to visit our website and factory to know more. We can supply the mould, R&D, custom and OEM/ODM service as per your request! Over 16 years Experiences in manufacturing & business, support customers with good service of rich experience which will save customer's working time to meet win-win!

Professional manufacturer



WELCOME TO ORDER



- Our aim is: "Top quality, Customer benefit firstly". Our company is based on "Credibility", aimed at "Win-Win" and took "service" as bounden duty. All products are strictly followed ISO9001:2008 and met CE, BSI, FCC, UL and RoHS standard.

## Production process



## Customer photoes



- Selectech Electronics is constantly improving itself at the spirit of building brand products and striding towards diversity and scale. Selectech has its own experienced technology engineers, modern management team, strong technical strength, advanced injection and testing machines and perfect quality control system to guarantee the quality, security and punctuality of delivery.



- Looking forward to cooperating with you!  
Selectech is your reliable and the best partner to start your business!

## Certificate



## Our advantages



## Why choose us



16 years OEM&ODM experience.



Custom-made products and packing as per customer's request.



Supply free cost samples.



All goods must be passed QA, QC, FQC inspect and test.



Quickly delivery time. Small qty have stock, Lots qty, it's about 20-25days.



We can supply drop shipping, cheapest and economic Express shipment.



We can collect all goods to ship for our customers.



Supply after-sales and quality guarantee service.



With FCC, CE, ETL, ROHS quality certification.



With ISO9001 Quality management system.



## Shipping choice



0000 00

1. 000 000 0 000?

0, 00 00 00 000 00000. 000 00000 0000.

2. 00 000 00000?

000 T/T, 000, 000 0000 0000000. 0000, 00 00 00.

3. 000 000 00 00 0 000?

000 00000 000 000 00 0 2400 000 000 000000. 000 00 000 00 00000 000000.

000 00 00 000 000 0 000 000 0000 0000000.

4. 000 00 0000 0 0 0000?

0. 000 USB 00 00 0 00 000 0 00 00 0000 000 000 000 00 00 0000 0000.

00. 000 000000 000000 000 000000 000 00 00000 00000 0 000 0000000. 000 00

000 000 000 000 000 00 0000. 00000 000, 00, 00000 000000 000 000 0000000.

0000 0000 0 0000. 000 00 000 000 000000000.

5. 00 000 00000?

000 EXW, FOB, CIF 00 000000. 00 000000 00 00000 00 000 0 0000.

6.000000 000000 000 000?

000 00 000 00000 0000 00 0000. 500pcs0 000000 0 8 0000000.

000000 000 000000 000000 200 000 000 00000 00 0000.

00.



7. 公司如何管理其供应链？

公司通过建立严格的供应商准入制度，定期对供应商进行考核，确保供应链的稳定性和产品质量。同时，公司还通过优化库存管理，降低物流成本，提高供应链效率。

8. 公司如何管理其研发项目？

公司采用OEM/ODM模式，与多家知名制造商合作，确保产品的质量和交货期。同时，公司还设立了专门的研发团队，负责产品的研发和创新，不断提升公司的核心竞争力。